



Specialized system for organic solvents



Alpha 3-4 LSCbasic



Product highlights

Special system Alpha 3-4 LSCbasic

4 kg

-105 °C

- ▶ **Coated real-glass chamber**
borosilicate glass for higher safety
- ▶ **Flow tube**
for optimum vapour flow
- ▶ **Deep-seated ice condenser**
to improve performance


Higher safety for the freeze drying of organic, solvent-containing samples

- Unique system set-up especially for organic solvents
- Flexible use for a wide range of different applications with organic solvents thanks to two fixed packages
- Ice condenser temperature of -105°C with natural refrigerants
- Particularly easy operation with the LSCbasic control system
- Capacitive vacuum sensor
- Controlled vacuum and solvent-resistant vacuum pump
- EPDM seal suitable for solvents
- Integrated valve for inertisation



Fixed configurations

For detailed information please visit our website:
martinchrist.de/alpha-3-4

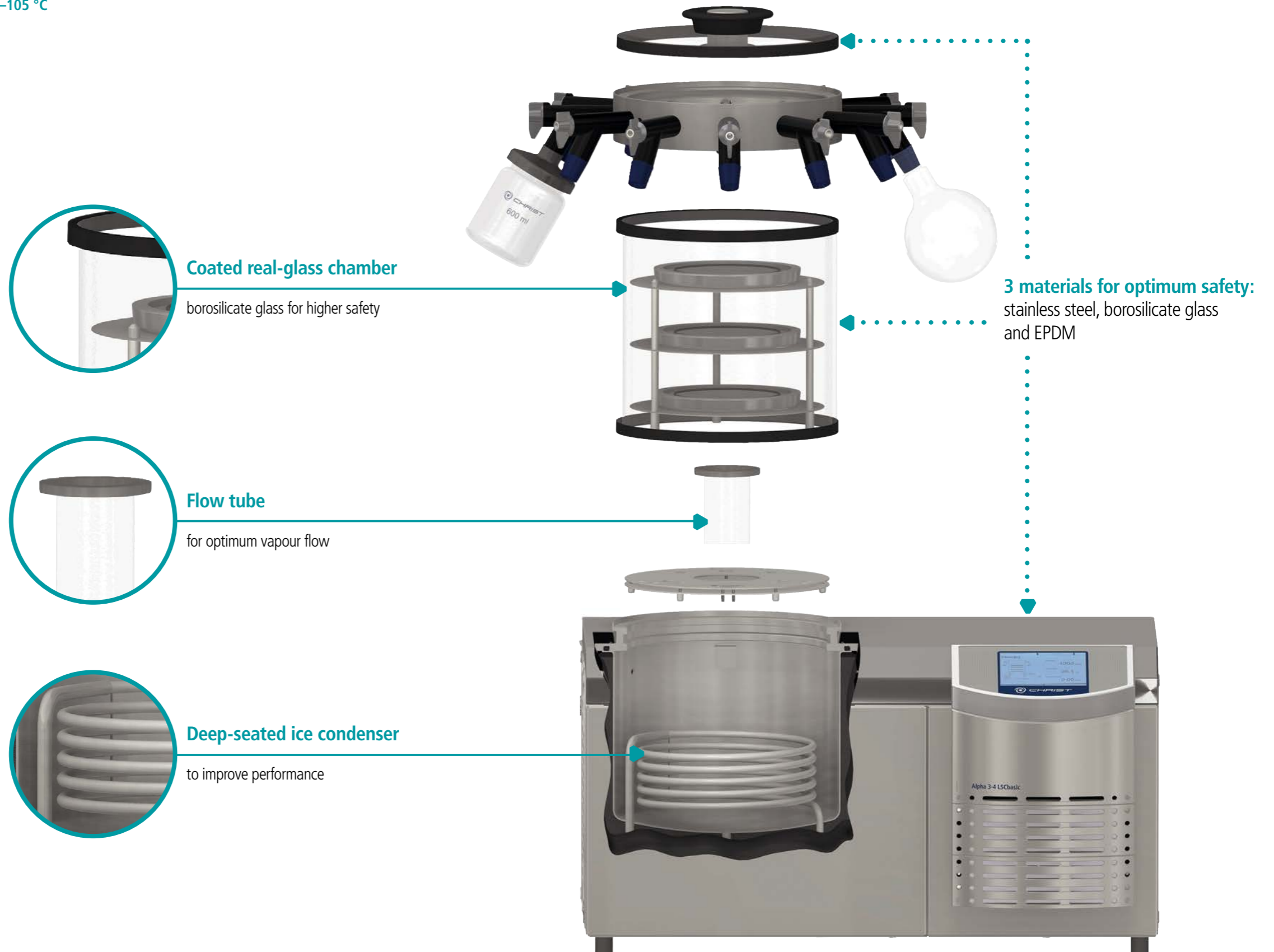



Martin Christ
Gefriertrocknungsanlagen GmbH

An der Unteren Söse 50
37520 Osterode am Harz, Germany

Tel. +49 (0) 55 22 50 07-0

Fax +49 (0) 55 22 50 07-12

info@martinchrist.de

www.martinchrist.de